Original Operating Instructions

S2

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dentalportal.info

This document is designed and released for the following group of persons:

- End users
- Authorized resellers
- Authorized service technicians

Distribution or duplication of this document is only allowed by written consent of vhf camfacture AG. This includes reproduction by presentation and broadcasting.

This document is an original document.

It is subject to technical changes and errors.

Contents

1 Legend	5
1.1 Explanation of safety message text	5
1.2 Classification of signal words	5
1.3 Symbols and hints	6
1.4 Instruction symbols	6
2 General safety messages	7
2.1 Supplemental directives	7
2.2 Safety messages	8
2.3 Intended use of the machine	10
2.4 Transport and storage	11
3 General information	12
3.1 About this manual	12
3.2 Machine components	13
3.3 Technical data	14
3.3.1 Basic system	14
3.3.2 Rotary axes	14
3.3.3 Spindle	14
3.3.4 Tool magazine and tool change	15
3.3.5 Blank changer	15
3.4 Manufacturing software	15
3.5 Scope of delivery	16
3.6 Sound emission	17
4 Installation	18
4.1 Placement requirements	18
4.2 Setting up the machine	19
4.3 Connecting the compressed air	20
4.3.1 Connecting the service unit	21
4.3.2 Setting the air pressure	22
4.4 Connecting the air extraction	22
4.5 Software installation	23

5	Running the machine	25
5.1	Operating the machine	25
5.2	Air extraction	25
5.3	Synchronous spindle SFS 300P	26
5.4	Tools	26
5.5	Blank changer	28
5.5.1	Equipping blank frames	28
5.5.2	Inserting blank frames	28
5.6	Procedure after a power failure	30
6	Maintenance and cleaning	34
6.1	Internal cleaning	34
6.2	Collet chuck cleaning	34
6.3	External cleaning (housing)	36
6.4	Main fuse change	37
6.5	Maintenance compressed air	37
6.6	Calibration	38
6.7	Definition of wear parts	39
6.8	Maintenance table	39
6.9	Inspection	40
7	Disposal	41
	Index	42
	EC Declaration of conformity	44

1 Legend

1.1 Explanation of safety message text

General warning text. Failure to comply with the information given in this text can lead to serious injury.

Nature and source of a hazardous situation

Possible consequences!

➤ Advice to avoid the hazardous situation

1.2 Classification of signal words

DANGER indicates a hazardous situation which, if not avoided, will result in death or serious injury.

WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury.

CAUTION indicates a hazardous situation which, if not avoided, could result in minor or moderate injury.

NOTICE is used to address practices not related to physical injury. It indicates situations that can lead to material damage of the product or the environment.

1.3 Symbols and hints

background information

- list

IMPORTANT

Operating instructions and other important information

HINT

Information to make work easier

1.4 Instruction symbols

★ requirement

M1. first action (M stands for manually)

S 2. second action in the software (S stands for Software)

M3. third action (M stands for manually)

► result

2 General safety messages

2.1 Supplemental directives

Read this manual before you start the machine. Follow the listed safety messages to avoid risks and possible major injuries.

The users have to be informed about the intended use of the machine. The safety messages indicated in this document have to be followed (see [chapter 2.3](#)).

Check the machine and especially the protective devices for possible damages. Damaged safety devices or parts thereof must, if not stated otherwise in the user manual, be repaired or replaced by authorised service personnel.

Besides a maintenance of the machine is executed, it may only be operated with a closed **working chamber door** and with the **safety interlock** being **activated**. It is **forbidden** to **disable** the present safety devices or to **set them out of use**.

Keep children and animals away from the machine. Don't leave the machine running unattended, except you comply with the requirements for it (see [chapter 2.3](#)). Use of the machine in any "service mode" may only be done by service personnel, that is trained by vhf to use the machine in the respective "service mode".

Please be aware of the information about the choice of tools and the adjustment of tool parameters (see [chapter 5.4](#)).

2.2 Safety messages

DANGER

Use of damaged cables

Danger to life due to an electric shock!

- Disconnect the machine from any electrical source and prevent the machine from being restarted.
- Contact the service department.
- Replace damaged cables only with original spare parts.

DANGER

Troubleshooting while the machine is running

Danger to life due to an electric shock!

- Consult the service department, before you clear the malfunction on your own.
- Disconnect the machine from any electrical source and prevent the machine from being restarted.

WARNING

Operation of the machine in service mode with open working chamber door

Risk of injury through cuts and bruises! Hazards through ejected chips!

- Operate the machine in “User” mode only unless you have been authorized by vhf to use other modes.
- Even if you are an authorized user, use the “Service Modes” only when necessary.
- While in any “Service Mode”: Do not reach into the working chamber while the axes are moving or during machining.
- While in any “Service Mode”: Everyone within reach of the machine must wear protective eye wear.

WARNING

Processing harmful material

Respiratory diseases due to inhalation of harmful materials!

- Process material only when the air extraction is running.
- Only use suited materials, that are harmless when the air extraction is working.
- Use a vacuum cleaner with an extra-fine particle filter.

 **WARNING**

Loud noise from the machine

Hearing damage due to loud noise!

- Change the milling conditions. Ensure that the blanks are fixed properly, the condition of the tool and check the material you use.
- Wear ear protection during processing, if loud noises can't be avoided.

 **WARNING**

Open pneumatic connections

Danger through loose pneumatic parts conducting compressed air!

- Separate the machine from the compressed air supply.
- Consult the service department.

 **CAUTION**

Disorder in the working environment

Trip hazard!

- Keep the working environment clean.
- Store the extra equipment like the vacuum cleaner or the wet grinding module in a safe place.

 **CAUTION**

Ergonomics in the working environment

Subsequent damage due to constant malposition!

- Adjust the working environment ergonomically.
- Pay attention to keep the seat height, the monitor position and the exposure to light ideally.

 **CAUTION**

Inserting milling tools

Cuts through sharp tool cutting edges!

- Grab milling tools only at the shank, not at the tool tip.
- Take care not to cut yourself or others while handling with the tool.

2.3 Intended use of the machine

The machines of the Impression S models have been designed for easy to medium milling works in the field of dental technology. It has to be considered that the machine cannot withstand all powers that can possibly occur during operation.

Only use the machine commercially.

The machines are designed to process the following materials:

- wax
- most of the plastic materials (e.g. PMMA)
- composites
- zirconium oxide
- non-precious alloys (CoCr-based)

The machine can get damaged, if it is used for another purpose than described above. The protective devices **may not** be **avoided** or **turned off**.

The instructions of this manual have to be followed and the maintenance procedures have to be executed as described.

The machine may only be operated with original equipment from vhf, to keep up the safety of the product and the warranty of the machine. The user carries the risk of damage to the machine if non-authorized equipment is used!

Wet grinding is exclusively permitted by using the wet grinding module from vhf.

The warranty expires, if you modify the machine or take off the housing without a written agreement of the vhf service department!

Unattended operation

The machine may run unattended, if the national and local laws are regarded. Besides that the requirements of the respective insurance company have to be met. The following requirements have to be met in any case, if the machine is run unattended:

- The working chamber of the machine has to be cleaned completely.
- Unauthorized persons may not have access to the machine.
- Don't change between dry processing and wet grinding.
- During wet grinding precautions have to be met, if fluid should leak.
- The room, in which the machine is standing, needs to have an automatic fire detection system.

Objects

Before creating jobs, verify if the objects being prepared may be utilized at the place of use according to local and/or national regulations or other authorized organisations or entities (e.g. professional associations, health authorities). In particular, verify if the material is approved for the machined object type and if the object type is designed in accordance with applicable regulations. Neither the manufacturing software nor the machine will inform you about possible regulatory infringements, but will execute jobs in accordance with the preferences and materials set by the user.

Verify that each object type and each material in your jobs are authorized manufacturing materials. If mandated by local or national regulations, obtain relevant authorization from the responsible organisation or entities (e.g. professional associations, health authorities).

Only manufacture objects which correspond to the object types that you can select in the import view in DentalCAM. While you can import/manufacture any other objects as well, neither the manufacturing software nor the machine are designed for these other objects and should not be used in this way.

Do not manufacture implants or parts of objects that are designed to have contact with implants. These parts include parts of two-part abutments which contain the connection geometry for the implant. Do not manipulate the connection geometry of prefabricated abutments ("prefab abutments") and you must always check finished objects for accurate connection geometries (i.e. that connection geometries of finished jobs have not been damaged).

2.4 Transport and storage

- The machine always has to be transported in an upright position.
- The machine has to be transported individually at all times. Several machines may not be piled up.
- For unpacking and positioning the machine two people are necessary. Use the shell type handle on the left and on the right side beneath the machine, to carry the machine.
- Ambient temperature (storage/transport): between -20°C and 60°C
- Relative air moisture: max. 80%, non-condensing.

3 General information

3.1 About this manual

This manual will enable you to use the dental milling machine S2 Impression as well as the accessory equipment safely and reliably. In this manual we want to advise you especially about the possible sources of danger that may come from the systems and to emphasise the correct use of the systems. Accordingly, we request you to read this manual attentively and to follow the advice given.

IMPORTANT

Keep this manual near the machine and ensure that everybody who is working with the machine has access to this manual.

We appreciate every feedback and suggestion for improvement, so that we can continuously improve our product and the associated documentation. If you have questions or suggestions, please refer directly to your reseller.

3.2 Machine components

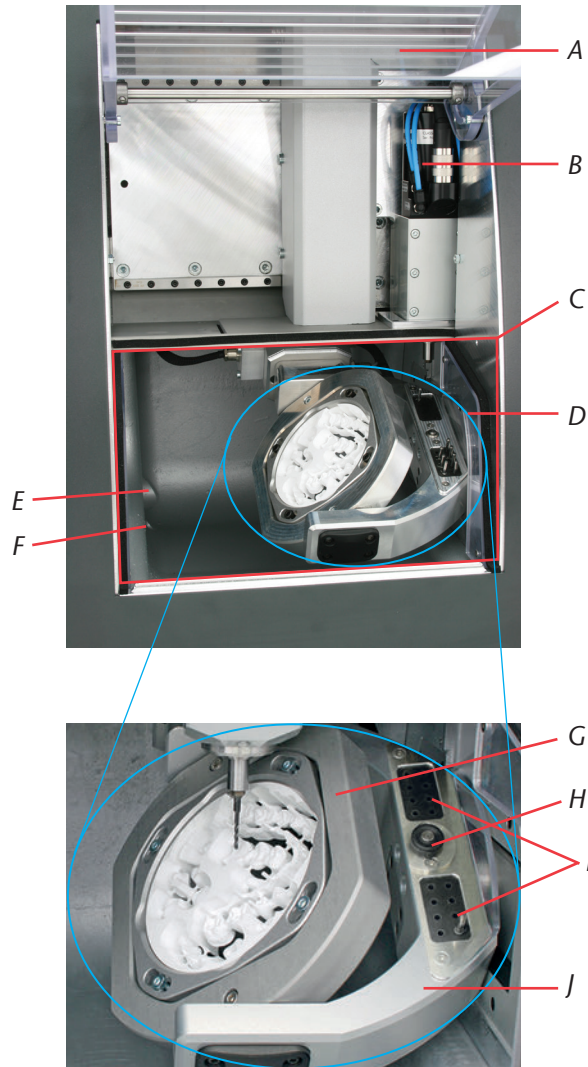
ILLUSTRATION 1: MACHINE COMPONENTS IMPRESSION S MODELS


- | | |
|---|---|
| A: working chamber door with safety contact and safety interlock | F: liquid drainage for wet grinding |
| B: synchronous spindle SFS 300P | G: A axis: rotation range of 360 degrees |
| C: working chamber | H: measuring key for tool length measurement |
| D: 4th/5th axis: blank holder for discs with \varnothing 98.5 mm and heights from 10 to 30 mm | I: tool magazine for 16 tools |
| E: opening for air extraction | J: B axis: rotation range of \pm 30 degrees |

3.3 Technical data

3.3.1 Basic system

- dimensions (W/D/H): approximately 692 x 445 x 540 mm
- five-axes mechanism; positioning range (x/y/z): 140 x 98 x 76 mm
- precise ball screw spindles for the 3 linear axes with 4 mm pitch
- repetition accuracy ± 0.003 mm
- complete housing of working chamber with working chamber door including a safety contact and safety interlock
- working chamber illumination
- prepared air extraction through hose connection and switching output at the back of the housing (requirements: external vacuum cleaner)
- underpressure sensor for monitoring the air extraction
- prepared wet grinding option with hose connection and switching output at the back of the housing (requirements: wet grinding module from vhf)
- disc dimensions: $\varnothing 98.5$ mm (± 0.2 mm), with step
- disc height: 10 mm to 30 mm
- weight S2 Impression: approximately 95 kg
- compressed air consumption of about 80 litres per minute

3.3.2 Rotary axes

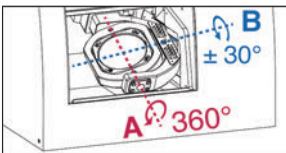- rotary axes free from backlash with Harmonic-Drive®
- rotation range: 360 degrees (A axis), ± 30 degrees (B axis)
- integrated blank holder for discs with $\varnothing 98.5$ mm

3.3.3 Spindle

- synchronous spindle with rotational speeds up to 60,000 RPM
- nominal power under continuous load (S1): 300 Watts
- maximal power (P_{max}): 600 Watts
- 4-fold bearing
- hybrid ceramic ball bearings
- pneumatic collet chuck 3 mm
- 3 liquid nozzles for wet grinding
- cone cleaning

3.3.4 Tool magazine and tool change

- automatic tool change
- tool magazine for 16 tools
- automatic tool length measurement and tool breakage monitoring via measuring key
- diamond-coated tools can be used
- compressed air monitoring for automatic tool change

3.3.5 Blank changer

- automatic blank changer with blank magazine, gripper and carriage to the working chamber
- holding up to 8 blank frames in the blank magazine
- safety contact at the blank magazine door
- pneumatically operated carriage cover between the working chamber and the blank changer

3.4 Manufacturing software

- i** To calculate the milling paths you need to create a file in the STL format from your CAD program. All important parameters for the respective material are included in the software package DentalCAM. Operation with external CAM programs is possible, if there is a post processor for the respective program. If you want to use an external CAM software contact vhf to clarify, if there is already a post processor for your CAM software or if a new post processor has to be written for your program. The machine is controlled with the included software DentalCNC.
- i** For operating the manufacturing software, please refer to the corresponding software manual.

3.5 Scope of delivery

ILLUSTRATION 2: SCOPE OF DELIVERY S2 IMPRESSION

- | | |
|--|--|
| A: 1 power supply cord | H: 1 key for emergency release of the working chamber door |
| B: 1 USB cable | I: 1 hose connection for external air extraction unit |
| C: 1 set of replacement screws for the blank holder and the tool measuring plate | J: 1 service unit for compressed air |
| D: 1 measuring pin | K: pneumatic hose \varnothing 6 mm |
| E: 1 drill 2.8 mm for the tool magazine inserts | L: 1 spindle service set |
| F: 2 tool magazine inserts | M: 1 micrometer |
| G: 1 hexagon socket screw key 2.5 mm to change the discs | N: 3 calibration discs |
| | O: 8 blank frames |

Further scope of delivery (not pictured):

- CNC milling machine S2 Impression
- manual S2 Impression
- software manual
- software license DentalCAM including control software. Please ask your reseller regularly, if new updates are available and update your software.

 If you are using an external CAM software only the control software DentalCNC is included in the scope of delivery.

3.6 Sound emission

 The sound emission varies heavily depending on the manufacturing material and the operating conditions. Change the milling conditions if the machine is too loud. Ensure that the blanks are fixed properly, the condition of the tool and check the material you use.

Wear ear protection during processing, if loud noises can't be avoided.

Measuring conditions:

- processed material: CoCr
- tool worn out
- measuring distance to sound source 1 m
- air extraction activated
- measurement according to ISO 3746, engineering method 3

Operating condition	Maximum A-weighted sound pressure level
processing (measuring conditions see above)	77 dB(A)
all other operating conditions (tool change, movement of the axes etc.)	<70 dB(A)

4 Installation

4.1 Placement requirements

The following requirements have to be considered to choose the right machine location:

- Firm and even surface, has to carry the weight of the machine (approximately 95 kg).
- Ambient temperature (storage / transport): between -20 °C and 60 °C
- Ambient temperature (operation): between 10 °C and 32 °C
- Relative air moisture: max. 80%, non-condensing.
- Machine location should be as dust-free as possible, since ambient air is sucked into the machine.
- Alternating current source with 100 – 240 V and 50 – 60 Hz
- Incoming air pressure maximum 8 bar
- Incoming liquid pressure of the wet grinding module maximum 3 bar
- The machine needs to have an adequate distance to the wall on the right side of the machine (seen from the front), so that the ventilation holes are not blocked (at least 10 cm distance). On the left side and on the back of the machine needs to be enough distance, so that you can reach all machine connections easily.

NOTICE

Plug the machine's power cord in a **separately fused** circuit current or ensure that no devices are connected that cause heavy voltage fluctuation when switched on. Heavy voltage fluctuations can disturb the control and cause a failure of the system.

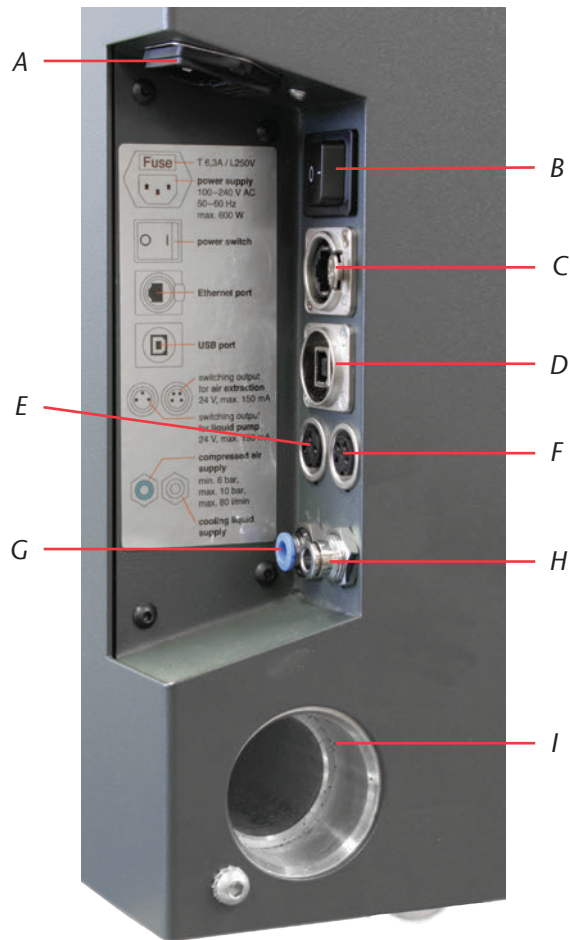
ILLUSTRATION 3: CONNECTION PANEL IMPRESSION S MODELS

- | | |
|---|--|
| A: power connection 100–240 V AC, 50/60 Hz including glass fuse T6,3A L250V | F: switching output for external air extraction unit |
| B: main power switch | G: pneumatic connection, Ø 6 mm, max pressure. 8 bar |
| C: optional Ethernet connection | H: liquid connection for external pump of the wet grinding module, Ø 6 mm, max. pressure 3 bar |
| D: USB connection | I: air extraction connection |
| E: switching output for external pump of the wet grinding module | |

4.2 Setting up the machine

IMPORTANT

Only use original cables from vhf to connect devices via the two switching outputs (see Illustration 3). If you are using a different USB cable than the one that is delivered, pay attention that it's not longer than 3 meters!

- M1.** Unwrap the milling machine. Please keep the packaging for possible further service shipments.
- M2.** Connect the service unit via the standard pneumatic connection with your compressed air supply (see also [chapter 4.3](#)).

IMPORTANT

The pneumatic hose (blue hose) must be inserted in the **blue** marked connection (on the left side at the bottom of the connection panel, see Illustration 3). It may **not** be inserted in the stainless steel connection, since the machine doesn't start without compressed air!

M3. Check the displayed air pressure value on the manometer. If the pressure is not between 7 and 8 bar, reset the pressure value (see [chapter 4.3.2](#)).

M4. Connect the CAM computer to the machine via the USB connection.

M5. Put the vacuum hose with the hose connection in the air extraction connection at the side of the machine (see also [chapter 4.4](#)).

M6. Connect the power cord to the machine.

IMPORTANT

M7. Plug the machine's power cord in a separately fused circuit current or ensure that no devices are connected that cause heavy voltage fluctuation when switched on (e.g. air compressor etc.).

M8. Close the working chamber door.

M9. Start the machine with the main power switch.

4.3 Connecting the compressed air

The spindle SFS 300P needs a pneumatic connection. The spindle uses the compressed air to control the collet chuck during the automatic tool change. It is also used to produce sealing air, that prevents dirt from entering into the spindle. The air consumption of the machine is 80 litres per minute.

The scope of delivery contains a service unit for the compressed air supply. You can use it to regulate the incoming compressed air to the machine. The built-in water separator containing a filter element prevents moisture or dirt from getting into the machine.

NOTICE

The incoming compressed air has to be **dry** and **oil-free** according to ISO 8573-1, since the water separator can filter only small amounts. Wet compressed air can lead to **bearing damage** and **electrical damage** of the spindle.

Connect the machine to the compressed air supply only via the provided service unit.

IMPORTANT

Air purity according to ISO 8573-1:

Solid particles	class 3	Filtration degree better than 5 microns for solid particles
Water content	class 4	Maximum pressure dew point +3 °C
Residual oil content	class 3	Maximum oil content 1 mg/m ³

The compressed air supply has to provide at least 6 bar during operation of the machine.

4.3.1 Connecting the service unit

The service unit can be mounted to the housing or connected directly to the compressed air supply. You need to make sure that the service unit is visible during operation so that a regular control is possible. In the following text the mounting to the housing is explained. If you connect the service unit directly to the compressed air supply, please regard only the steps **M3.** to **M5.**

NOTICE

The service unit **always** has to be installed **vertically** (see Illustration 4 on page 22). The water separator doesn't work, if the service unit is not installed vertically!

- M1.** Remove the two blind screws under the ventilation holes at the side of the machine.
- M2.** Mount the service unit vertically with the included cylinder head screws (see Illustration 5 on page 24).

 **CAUTION**

Connect the maintenance unit to the machine **before** you connect the maintenance unit to the compressed air supply. Otherwise risk of injury occurs through leaking compressed air or lashing pneumatic hoses.

- M3.** Cut a piece of the delivered pneumatic hose (Ø 6 mm) in the desired length.
- M4.** Connect the pneumatic connection on the right side of the maintenance unit to the pneumatic connection of the machine.
- M5.** Connect the compressed air supply to the standard pneumatic connection on the left side of the maintenance unit.

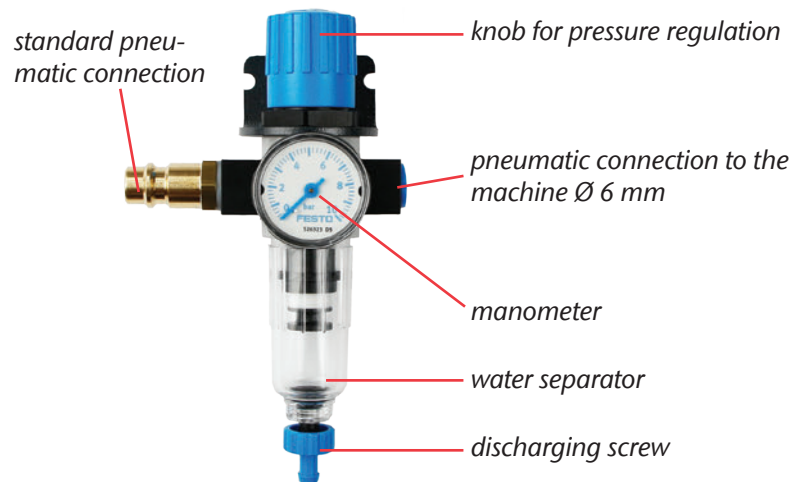
ILLUSTRATION 4: SERVICE UNIT: AIR PRESSURE SETTING AND CONTROL

4.3.2 Setting the air pressure

HINT

Setting the air pressure is only necessary, if the displayed pressure on the manometer is not between 7 and 8 bar.

★ The maintenance unit is connected to your compressed air supply.

M1. Pull the knob on top of the service unit a little bit.

M2. Set the pressure (turning the knob toward “+” increases the pressure, turning the knob toward “-” decreases the pressure), until the displayed value is at least 7 bar.

M3. Push the knob down again.

► The air pressure is set and cannot be changed inadvertently.

4.4 Connecting the air extraction

IMPORTANT

The vacuum cleaner you use has to be suitable for your application (**M-class filter** for processing harmful material etc.).

- Empty the dust bag of your air extraction regularly and replace the possibly existing filter in regular intervals.
- Read the manual of your air extraction carefully and follow the advice given.

M1. Plug the hose connection with the air extraction hose in the air extraction connection at the side of the machine.

M2. Connect the optional switching unit to the switching output at the side of the machine (see Illustration 3 on page 19).

M3. Connect the air extraction to the power supply respectively connect it to the switching unit and connect the switching unit to the power supply (see Illustration 5 on page 24).

4.5 Software installation

★ The CAM computer is running and connected to the machine.

S 1. Start the installation file of the software. You can receive the installation file from your vhf reseller.

S 2. Unzip the installation file.

S 3. Start the file setup.exe and follow the advice of the installation program.

Further advice for the software installation is given in the software manual.

Please ask your reseller regularly, if new updates are available and update your software.

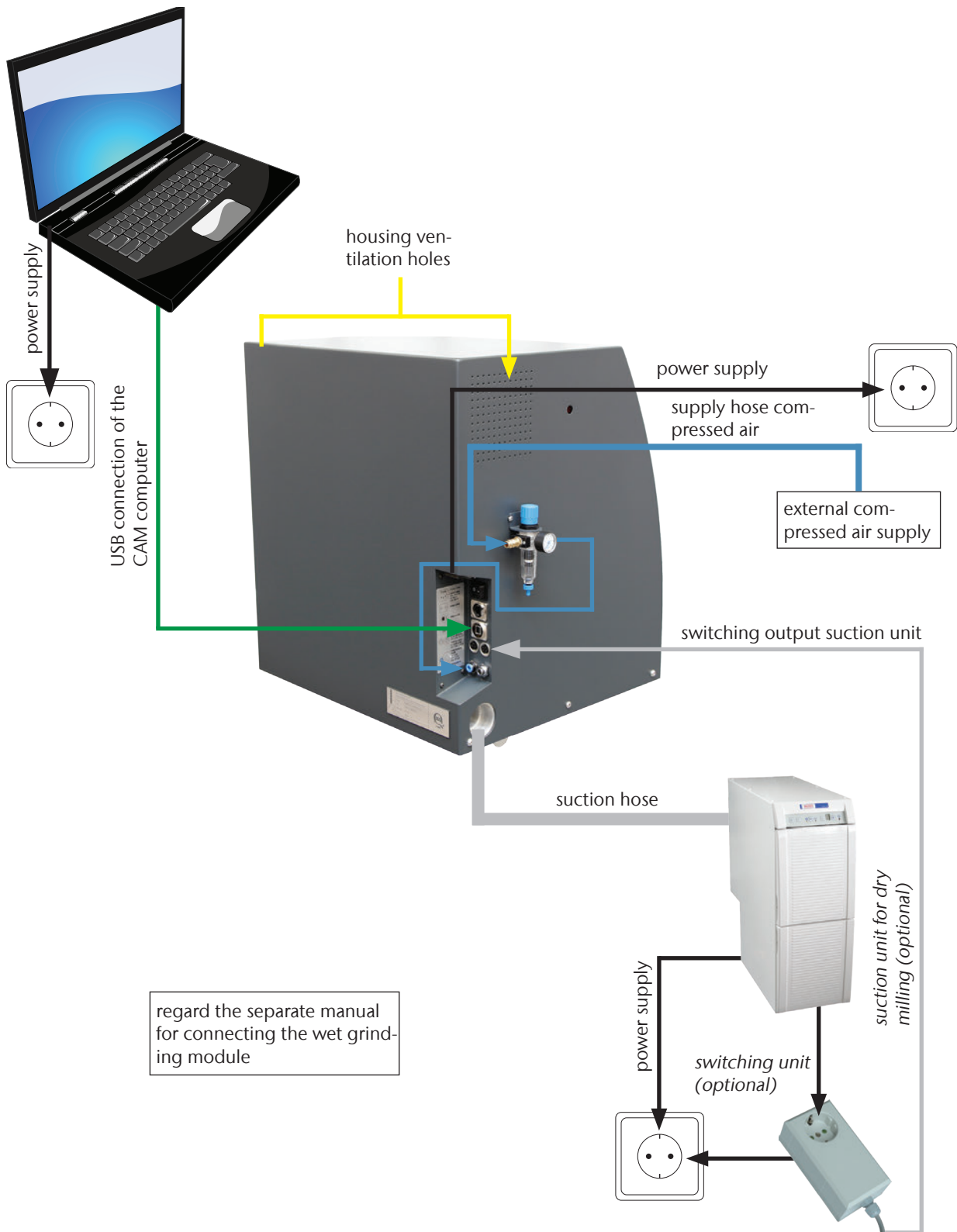
ILLUSTRATION 5: OVERVIEW OF THE MACHINE CONNECTIONS IMPRESSION S MODELS

5 Running the machine

5.1 Operating the machine

Make sure that all cables and hoses have been connected correctly and turn on the machine with the main power switch on the side of the machine.

For a milling job you only need the blanks which you want to use as well as the appropriate tools. As the machine is operated through the software, you can find further details on handling the machine in the software manual.

The machine is equipped with a safety contact and a safety interlock for the working chamber door. When the safety interlock is activated, the working chamber door cannot be opened during a milling operation.

IMPORTANT

The machine may **only** be operated when the **safety interlock** is **activated**.


5.2 Air extraction

NOTICE

Operate the machine only with the **air extraction turned on** and make sure, that the air extraction opening is not blocked. If the air extraction is not working properly, material chips can reach the sensitive parts of the machine like the bearings or the ball screws and damage them!

The models of the Impression S models are prepared for air extraction with an external industrial vacuum cleaner. For this purpose, the machine is equipped with an opening for the hose connection as well as a switching output for automatically switching the suction unit at the housing side.

For operating the air extraction unit you need the provided hose connection and an external vacuum cleaner.

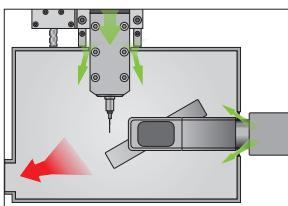The machine is equipped with an underpressure sensor for monitoring the air extraction. Make sure that the air extraction is running and the air extraction connection is not blocked. The integrated **working chamber sealing air concept** prevents milling dust and chippings from entering into the mechanics and electronics of the machine and reduces wear and maintenance efforts.

5.3 Synchronous spindle SFS 300P

WARNING

For the operation, installation and maintenance of the spindle follow the regulations for safety at work. Improper handling or operations differing from the intended use of the machine greatly reduce the safety of usage!

NOTICE

The synchronous spindle SFS 300P is a high precision device that may never be subject to any form of raw force like strikes!

Keep in mind that the values for the immersion depth and the feed rate are not too high during operation. This can lead to damage of the blanks and tilting of the tool. By that the spindle can take irreparable damages. If you work with the DentalCAM software, the software automatically calculates the optimal parameters.

Do not use unbalanced tools (single tooth cutters, gravers) at high rotational speeds. Such an imbalance puts a great strain on the ball bearings of the spindle, by that the bearings can get damaged.

5.4 Tools

CAUTION

| **Never** touch a tool at the tip. Otherwise cutting injuries may occur.

HINT

Usage of original tools is recommended, since the tools are designed by vhf especially for the machine and the designated milling operation.

IMPORTANT

The tools used for cobalt-chrome blanks are shorter than tools for other materials. Therefore only cobalt-chrome blanks with a **thickness of maximal 18 mm** can be used on the machine!

If you should use tools from other manufacturers, please regard the following specifications.

Only tools with a shank diameter of 3 mm and a total length of up to 40 mm may be used. For storing the tool in the tool magazine an unmachined shank with a length of at least 7 mm is required. The clamping length in the collet chuck should be between 13 and 14 mm.

NOTICE

The tool needs a sufficiently large chamfer at the tool shank, since the collet chuck can get damaged otherwise.

A retaining ring has to be installed as a stop ring into the existing slot. Install only retaining rings according to DIN 471-A3!

Make sure that the tools are put straight into the tool magazine. All inserted tools must suit the **positions stored in the software**. Otherwise the **tool** or the **blanks** may get **damaged**.

Don't use tools, where the cutting edge diameter is higher than the shank diameter (\varnothing 3 mm). **Don't** use tools with a cutting edge diameter higher than 2.5 mm for the tool magazine, otherwise the tool magazine inserts get damaged.

5.5 Blank changer

5.5.1 Equipping blank frames

- M1.** Loosen the 4 screws of the blank frame (see Illustration 6) with the delivered hexagonal socket screw (2,5 mm).
- M2.** Insert a disc.
- M3.** Apply the fixing bows.
- M4.** Tighten the screws hand-tight.

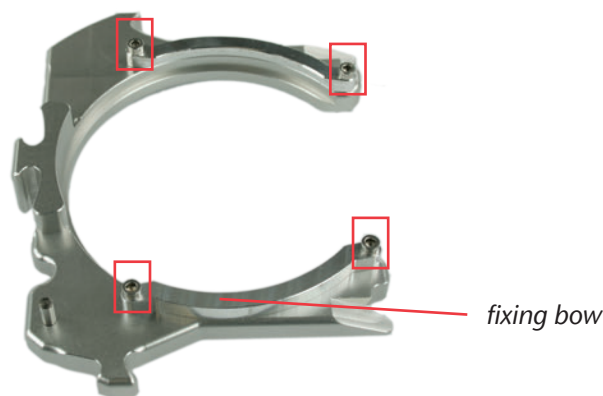
ILLUSTRATION 6: BLANK FRAME S2 IMPRESSION

5.5.2 Inserting blank frames

Further advice for operating the blank changer is given in the software manual.

IMPORTANT

The blank magazine door is equipped with a safety contact. The blank changer can only move when the blank magazine door is closed.

★ The machine is turned on, the working chamber door and the blank magazine door are closed.

S 1. Choose the desired blank magazine slot in the software.

► The blank magazine turns until the respective blank frame is at the change position.

M2. Open the blank magazine door.

M3. Grab the blank frame in the grasp area across from the clamping bolt (see Illustration 7 on page 29 in the lower right corner).

CAUTION

Crushing hazard through movement of the clamping bolt! Grab the blank frame only in the grasp area. **Don't** grab the blank frame in the area near the clamping bolt!

M4. Push the frame release and keep it engaged.

► The clamping lever opens to release the frame at the change position.

HINT

You can push the frame release and insert the blank frame with one hand only.

M5. Insert the blank frame at the change position (vertical position).

M6. Disengage the frame release.

- The clamping lever closes and fixates the blank frame.

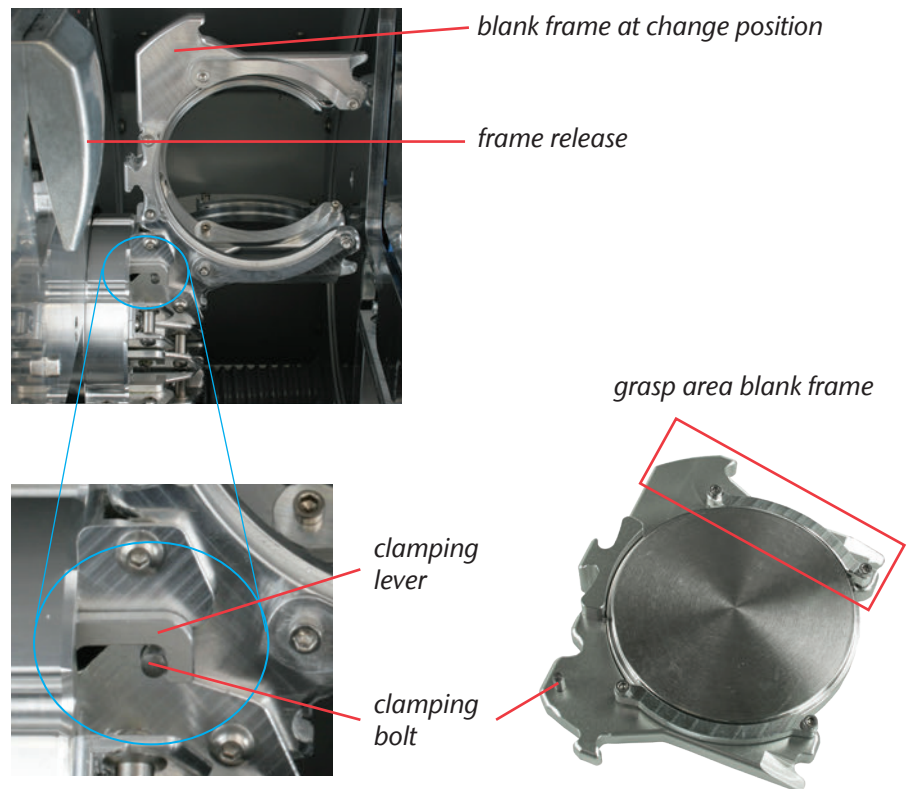
ILLUSTRATION 7: EQUIPPING THE BLANK MAGAZINE

5.6 Procedure after a power failure

If there is a short power failure while the blank changer is not in use, it is enough to turn the machine back on again. If a power failure occurred during processing, follow the advice of the software, so that the machine can reference.

If there is a longer power failure or a malfunction of the power wire you can deactivate the safety interlock with an emergency key to take the blank out of the machine.

IMPORTANT

If the power supply is interrupted while the blank changer operates, it's not enough just to turn the machine back on again. In this case the blank changer produces an emergency stop signal, before the machine references. This is a safety feature to protect the machine from getting damaged.

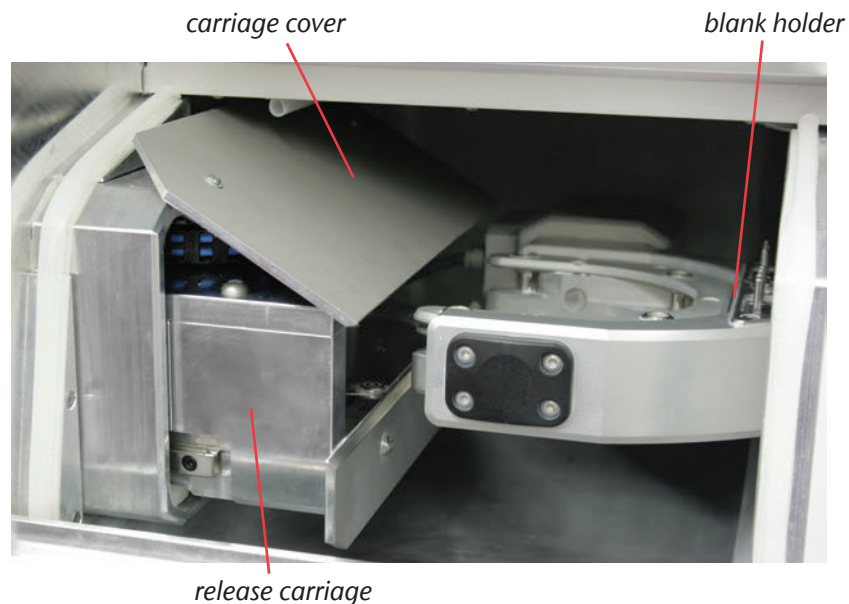
ILLUSTRATION 8: RELEASE CARRIAGE S2 IMPRESSION

★ The machine had a power failure, while the blank changer was in operation.

M1. Turn off the machine at the main power switch and disconnect the machine from any electrical source and the compressed air supply.

M2. Open the blank magazine door.

- M3.** Lead the key for the emergency release of the working chamber door through the opening to the right of the blank changer LED (see Illustration 9).

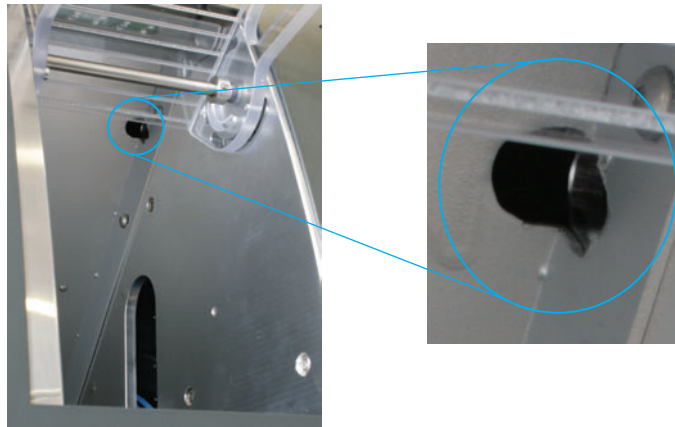
ILLUSTRATION 9: OPENING OF THE EMERGENCY RELEASE

- M4.** Deactivate the safety interlock of the working chamber door by turning the key counter-clockwise until the end position.

► The working chamber door can be opened.

- M5.** Open the working chamber door.

- M6.** Open the gripper (see Illustration 11 on page 33) and take out the blank frame, if it is not locked in the blank holder.

- M7.** Close the gripper.

HINT

If you don't reach the gripper, push the release carriage to the blank magazine until you reach the gripper.

NOTICE

Push the gripper only against the left stop when the gripper is closed to protect the machine from getting damaged (see Illustration 11 on page 33)!

- M8.** Hold the cover plate up with one hand and push the release carriage to the blank magazine until the end position with the other hand (see Illustration 10).

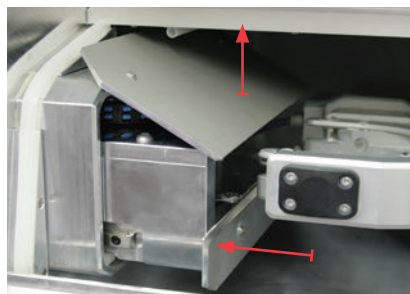
ILLUSTRATION 10: RELEASE CARRIAGE IN WORKING CHAMBER

M9. Push the gripper away from the working chamber until the end position at the side panel.

M10. Open the gripper.

- The release carriage is in the basic position.

M11. Close the working chamber door.

M12. Activate the safety interlock of the working chamber again by turning the key clockwise until the end position.

IMPORTANT

Make sure that the safety interlock is activated again. The machine may **only** be operated when the safety interlock is **activated**.

M13. Turn the machine back on again.

- The machine references and is ready for use again.

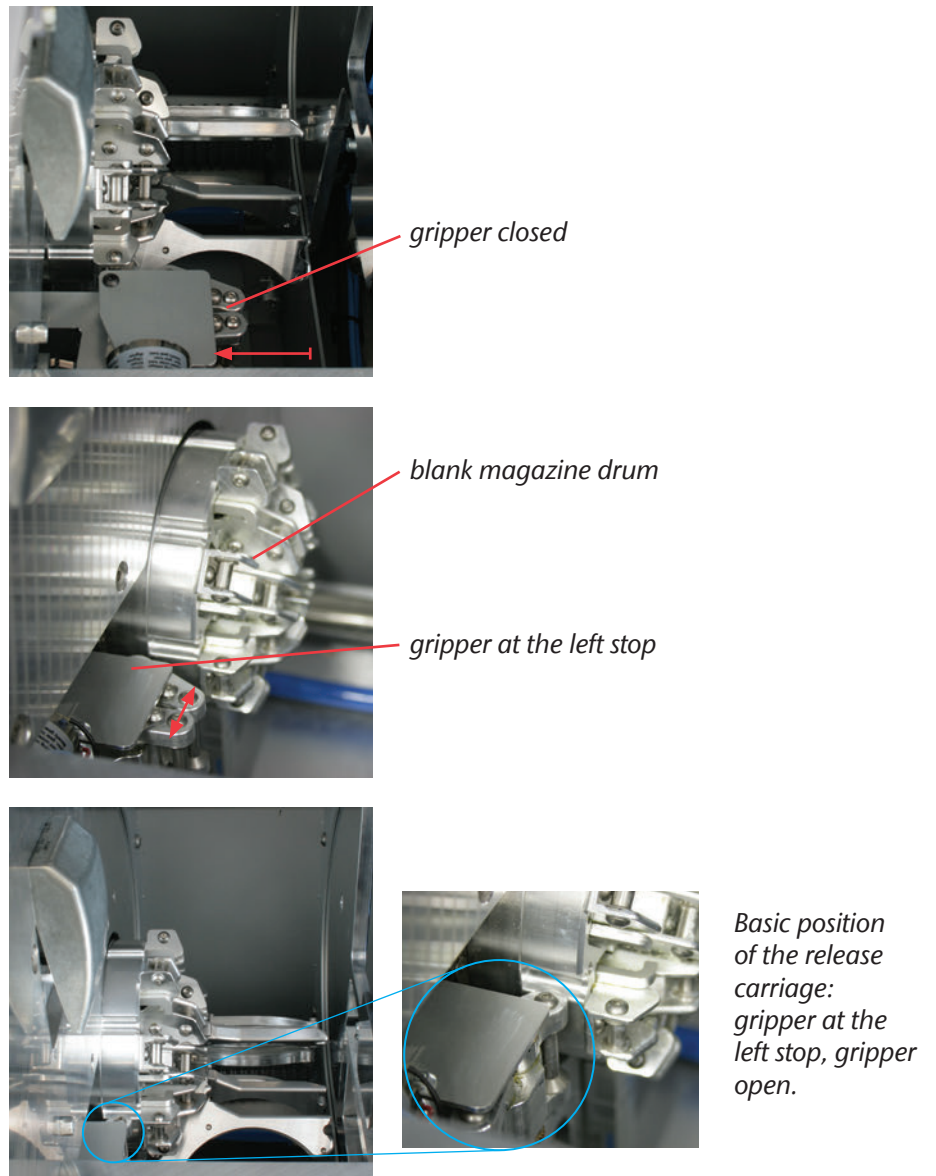
ILLUSTRATION 11: BASIC POSITION RELEASE CARRIAGE

6 Maintenance and cleaning

6.1 Internal cleaning

NOTICE

Never use **compressed air** for internal cleaning since material chips can reach the sensitive components like the linear guides or the spindle bearings.

Clean the machine regularly to prevent material chips from getting to the linear guides and other sensitive machine parts. Use only a vacuum cleaner or a brush for internal cleaning.

Clean the fixing device and the measuring key as well with a vacuum cleaner and a brush.

HINT

If you only use a brush for cleaning, turn on the air extraction in the software, so that the chips are removed instantly.

To clean the blank changer remove all blank frames and clean them. Clean the interior of the blank changer with a vacuum cleaner and a wet cloth. The best way to reach the release carriage is by lifting the carriage cover in the working chamber.

6.2 Collet chuck cleaning

HINWEIS

Damaging of the spindle when cleaning with compressed air

If you clean the collet chuck with compressed air, the spindle bearings can get damaged.

- Clean the collet chuck exclusively with the appropriate service set.

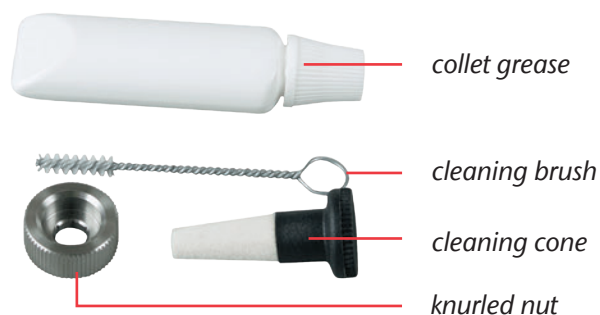
ILLUSTRATION 12: SCOPE OF DELIVERY SPINDLE SERVICE SET

- ★ You have the delivered spindle service set ready.
The working chamber door of the machine is closed.
- M1.** Open the collet chuck with the “Release tool” function of DentalCNC.
- M2.** Open the working chamber door and remove the tool from the spindle if necessary.
- M3.** Put the knurled nut to the spindle with one hand. With the other hand insert the measuring pin through the knurled nut into the collet chuck and keep holding it.
- M4.** Unscrew the collet chuck with the knurled nut.
- M5.** Remove the measuring pin from the collet chuck. Then put it aside within reach together with the knurled nut.
- M6.** Clean the inner cone of the spindle with the cleaning cone of the service set.
- M7.** Clean the collet chuck with the brush of the service set. Make sure that the 3 buffers in the collet chuck do not fall out during cleaning.

HINWEIS

Damaging of the spindle when using the wrong or wrongly applied grease

If you use grease that is unsuitable for the spindle or if grease gets into the longitudinal slots of the collet chuck, the spindle can get damaged.

- Make sure that no grease gets into the longitudinal slots of the collet chuck.
 - Only use a very small, about pinhead-size amount of the grease.
 - Only use the provided collet grease of the service set.
- M8.** Put a small amount of the collet grease on the index finger and smear it with the thumb.
- M9.** Apply the smeared collet grease to the flanks of the collet chuck.

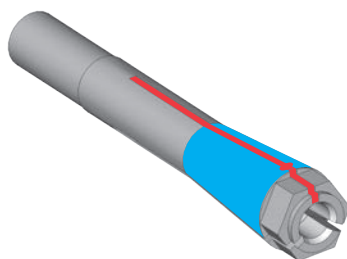
ILLUSTRATION 13:

GREASING THE COLLET CHUCK; SURFACE TO GREASE MARKED BLUE; SLOT THAT MAY NOT BE GREASED MARKED RED

- M10.** Insert the measuring pin into the collet chuck with one hand and keep holding it. Screw the collet chuck with the knurled nut in your other hand tightly into the spindle.

IMPORTANT

In any case, turn the knurled nut as far as it goes so that the collet chuck sits tightly in the spindle. Otherwise rotational imperfections can occur during operation which worsens your processing results.

M11. Remove the measuring pin from the collet chuck and store it together with the other components of the spindle service set.

M12. Close the working chamber door.

6.3 External cleaning (housing)

The cleaning of powder-coated surfaces should generally be done with a soft dry cloth first. If pollution cannot be removed this way, the cloth can be moistened and a pH-adjusted cleaner can be used if necessary.


IMPORTANT

Basically any contact between the powder coatings and alkaline or acid substances has to be avoided. Especially metallic powders show a high sensitive reaction. Please take special care of the affixed machine labels, so that you don't remove them while cleaning the housing. They are highly sensitive to rubbing and strong detergents.

HINT

If the usage of a special cleaner is necessary to eliminate certain pollutions we recommend to check the applicability of the cleaner at a hidden place of the housing first.

6.4 Main fuse change

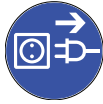★ The main fuse of the machine is damaged.

- M1.** Turn off the machine at the main power switch and disconnect the machine from any electrical source and the compressed air supply.
- M2.** Remove the cable of the power supply at the connection panel.
- M3.** Remove the cover of the fuse (bottom part of the power supply socket, see Illustration 14)
- M4.** Remove the damaged fuse and replace the damaged fuse by a new fuse.

IMPORTANT

Only use a glass fuse of the type T6.3A L250V as replacement fuse. A new main fuse is available as spare part at vhf.

- M5.** If you don't have a replacement fuse ready, take the replacement fuse from the right side of the fuse cover and put it into the left side of the fuse cover.
- M6.** Remount the fuse cover.

ILLUSTRATION 14: MAIN FUSE CHANGE

6.5 Maintenance compressed air

Please check regularly if there is liquid in the water separator of the service unit and empty it if necessary by opening the discharging screw (see Illustration 4 on page 22). The condensate will be blown out under pressure. Condensate in the separator usually points to compressed air being not dry enough. Please note our guidelines for air purity in [chapter 4.3](#).

IMPORTANT

In case of strong dirt accumulation the filter element in the water separator has to be cleaned or exchanged. A strongly polluted filter element may cause a loss of air pressure.

- M1.** Separate the machine from the compressed air source.
- M2.** Unscrew the bowl of the water separator.
- M3.** Turn out the filter screw under the filter element.
- M4.** Pull out the filter element.

A new filter element is available as spare part at vhf.

M5. Insert the new or cleaned filter element and reassemble the water separator.

- i** If the compressed air complies with the standard ISO 8573-1 (see [chapter 4.3](#)), the filter element usually doesn't have to be exchanged at all. If the filter element is polluted, please control the air purity of your compressed air.

6.6 Calibration

The machine is already calibrated on delivery. As long as your milling results are correct, there is no need for a further calibration. If your results should become inaccurate, a calibration could become necessary. In such a case, please contact a service technician before recalibrating the machine.

A calibration is time consuming and can degrade your milling results in the worst case, if it is not done correctly. Please try to change the milling conditions first. Check the fixation of the blanks, the condition of the tool and check the material you use, before you calibrate the machine.

- i** Further advice for the software installation is given in the software manual.

6.7 Definition of wear parts

The machine and the extra equipment are warranted for a period of 24 months or 2000 operating hours, whatever comes first. The warranty covers defects in materials or fabrication as long as the regulations for using the machine in all documents are followed.

Of course, the warranty also covers wear parts as long as their failure cannot be attributed to the function-related abrasion. The wear parts that are mentioned below can already wear down within the warranty period due to their normal function. The average useful life of the wear parts can be seen in the following table. Use these values to determine operating costs, to plan your spare part stock as well as to create individual maintenance and service plans.

6.8 Maintenance table

Task	Recommended interval	Procedure/Utensils	Illustration spare part
Cleaning the working chamber	Once per day, if necessary more often	Cleaning with vacuum cleaner and brush (chapter 6.1)	—
Collet chuck cleaning	Once a week, if necessary more often	Spindle service set (chapter 6.2)	
Cleaning the housing	As necessary	Cloth and possibly detergents (chapter 6.3)	—
Calibration of the machine	Only if necessary (chapter 6.6)	CNC Software (see software manual)	—
Changing the filter cartridge of the service unit	With visible pollution or every 2 years*	Check compressed air supply (chapter 6.5)	
Inspection by service technician	After 2000 operating hours* or every 2 years* (chapter 6.9)	—	—

* These are recommendations as a guidance. Depending on the processing material and how well the machine is cleaned these values can differ.

Wear parts

Exchanging wear parts	Average useful life	Procedure/ Utensils	Illustration spare part
Exchanging the spindle bearings (requires spindle exchange by service technician)	After 2000 operating hours* or every 2 years*	Performed by service technician	
Exchanging the tool magazine inserts	After 500 operating hours*	CNC Software (see software manual)	
Exchanging the collet chuck	After 1000 operating hours*	—	

* These are recommendations as a guidance. Depending on the processing material and how well the machine is cleaned, these values can differ.

6.9 Inspection

vhf recommends to have a service technician perform an inspection regularly. Additionally, have an inspection carried out whenever the machine is being repaired.

7 Disposal

Disposal of the machine is taken over cost-free by vhf. The customer/user bears the costs for disassembly, transport and packaging of the machine. Contact vhf in advance of a final disposal to go through the further procedure.

Please regard the national and the local laws of the disposal location, if you dispose the machine on your own.

Index

A

Air extraction.....	25
Air extraction, connection	22
Air pressure, setting.....	22
Air purity	21
Ambient temperature	11, 18

B

Blank changer.....	28
Blank frame	28
Blank magazine slot.....	28

C

Calibration	38
Carriage cover.....	30
Change position (blank frame).....	29
Compressed air, connection.....	20
Compressed air, consumption.....	14
Compressed air, maintenance	37
Connection panel.....	19

D

Discs	
Dimensions	14
Mounting	28
Disposal	41

E

External cleaning	36
-------------------------	----

F

Feedback.....	12
Filter element	37
Frame release	29

I

Inspection	40
Intended use of the machine	10
Internal cleaning	34
ISO 8573-1	21

L

Legend.....	5
-------------	---

M

Machine components.....	13
Machine connections	24
Machine location.....	18
Main fuse.....	37
Main power switch	19
Maintenance table.....	39
Manual	12
Manufacturing software	15
Materials	10

P

Packaging.....	19
Placement conditions	18
Power failure	30

R

Release carriage	30
Rotary axes	14

S

Safety messages.....	8
Scope of delivery.....	16
Service unit, connection.....	21
Setting up the machine	19
Software installation	23
Software manual	15, 23, 38
Sound emission	17
Spindle.....	14, 26
Spindle service set.....	34
Storage.....	11
Supplemental directives	7

T

Technical data	14
Tool magazine	15
Tools	26
Transport.....	11

U

Underpressure sensor.....	25
---------------------------	----

W

Warranty.....	10, 39
Wet grinding module	14
Working chamber.....	13
Working chamber sealing air concept.....	25

EC Declaration of conformity

according to EC directive for machinery 2006/42/EC Annex II A

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hereby declare expressly, that the

Machine: CNC milling machine
Type: CAM 5-S1 Impression ; CAM 5-S2 Impression
Serial No.: S1ID100000000 – S1ID199999999 ;
S2ID100000000 – S2ID199999999

fulfils all the relevant provisions of the following directives:

- 2006/42/EC Machinery directive
- 2014/30/EU EMC directive

References of the applied harmonised standards according to article 7 passage 2:

- EN 614-1:2006+A1:2009
- EN ISO 12100:2010
- EN 13128:2001+A2:2009
- EN ISO 13849-1:2015
- EN ISO 13849-2:2012
- EN 60204-1:2006
- EN 61326-1:2013
- EN 61326-2-1:2013
- EN 61000-3-2:2014
- EN 61000-3-3:2013

References of other standards:

- IEC 61010-1:2010 +A1:2016

The manufacturer undertakes to transmit, in response to a reasoned request by the national authorities, relevant information on the machinery electronically. Person being established in the community, who is authorised to compile the technical file:

Dipl.-Ing. (FH) Frank Benzinger
Vorstandsvorsitzender / Chief Executive Officer (CEO)
vhf camfacture AG
Lettenstraße 10
D-72119 Ammerbuch

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